



**For safe and
uninterrupted
laparoscopy**

Essential Laparoscopic Suite





As a large-scale provider of operating theatre-related products, laparoscopy is a core part of Medline's expertise. Our suite of fundamental laparoscopic products enables you to perform minimally invasive procedures, with flexibility and consistency. It is a comprehensive solution that facilitates the achievement of continuity in the operating theatre while focusing on safe surgery.


MED-LAP ORIS Insufflation Tubing Set

For proper CO₂ insufflation

The ORIS Insufflation Tubing Set completes Medline's laparoscopic portfolio. By filtering and insufflating clean CO₂, it allows clinicians to effectively establish pneumoperitoneum prior to and during laparoscopic procedures.

Medline offers a line of items that are fully compatible for use in conjunction with different insufflators.

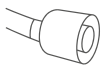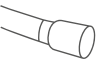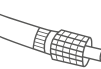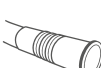Image	ORIS Insufflation Tubing Set	Compatibility range
	ORISISO with ISO 15 M/22 F connector	Karl Storz Richard Wolf Stryker ConMed
	ORIS15 with 15-mm male connector	Karl Storz Richard Wolf
	ORISLL with Leur-lock connector	Olympus
	ORISPUSH with flexible female connector	Stryker

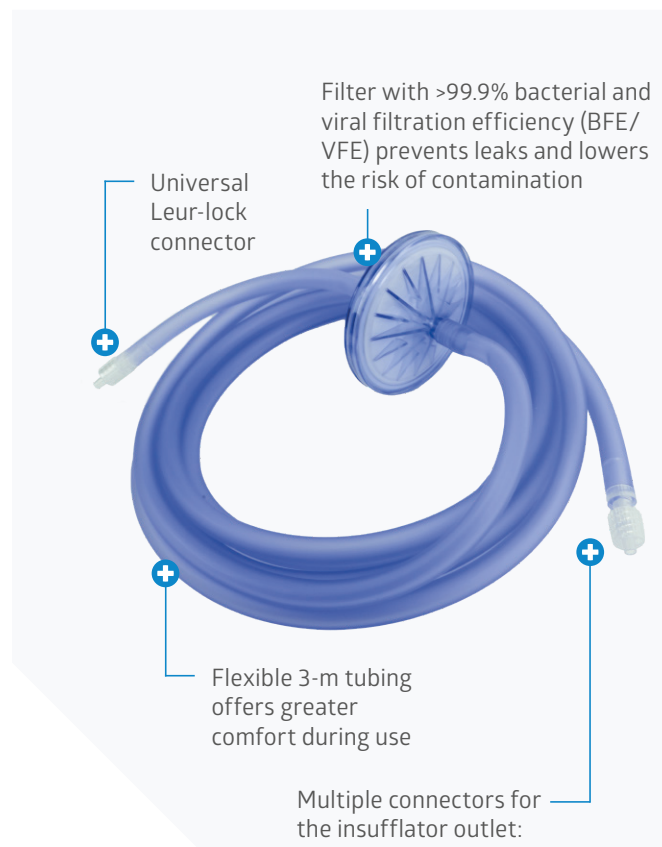

Benefits and features:

- Offers flexible tubing and a highly secure filter
- Filters with 99.9% BFE and VFE
- Fully compatible range



ORISLL
Luer-lock connector



ORISPUSH
Female connector




ORISISO
ISO 15 male/22 female connector

Item code	Description	Length	Packaging
ORISISO	Operating theatre insufflation set with ISO 15 M/22 F connector	3 m	25 ea/cs
ORISS15	Operating insufflation set with 15-mm connector	3 m	25 ea/cs
ORISLL	Operating theatre insufflation set with Luer-lock connector	3 m	25 ea/cs
ORISPUSH	Operating theatre insufflation set with female connector	3 m	25 ea/cs

MED-LAP Veress Needle

For safe establishment of the pneumoperitoneum

The insufflation needle is equipped with a retracting, spring-loaded blunt stylet and is designed for simple and safe establishment of pneumoperitoneum during laparoscopic procedures.

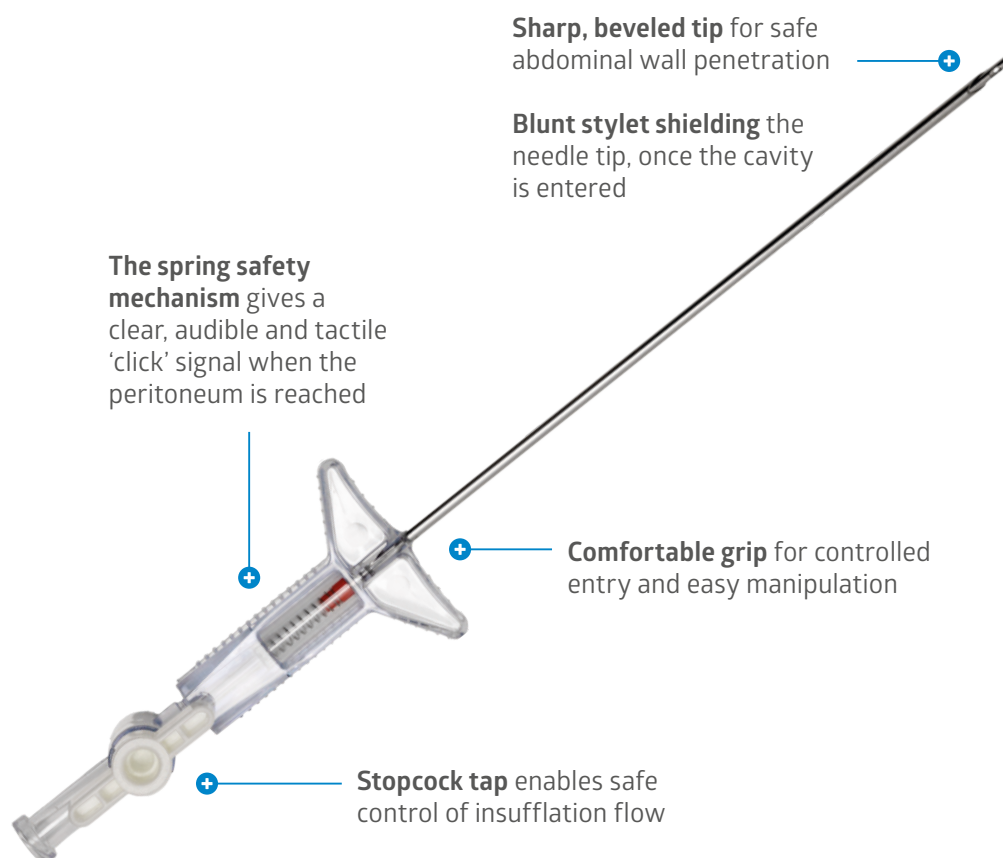

Benefits and features:

- Two working length options: 120 mm and 150 mm
- Standard luer lock connector to connect insufflation tubing
- Standard 14-G diameter

Item code	Description	Length	Gauge	Insufflation eyes	Packaging
ORVN12	Veress insufflation needle	120 mm	14	2	12 ea/cs
ORVN15	Veress insufflation needle	150 mm	14	2	12 ea/cs

MED-LAP Anti-Fog Solution

For preventing foggy laparoscopic lenses

The Anti-Fog Solution prevents laparoscopic lenses from misting, thereby providing an unobstructed view during laparoscopic procedures.



+ 6-ml volume that lasts for the entire length of a procedure



+ Radiopaque and adhesive sponge, for convenient placement in operating field



Benefits and features:


- Non-abrasive lint-free applicator sponge
- Ready-to-use kit: bottle and sponge
- Screw-cap bottle with mild alcohol solution

Item code	Description	Packaging
NONFB100	6 ml	20 ea/cs

MED-LAP Suction-irrigation Set

For keeping the working field clean and clear of liquids and smoke

This suction-irrigation system enables easy and efficient aspiration or irrigation, in order to keep the abdominal cavity clear of blood, liquids and smoke during laparoscopic procedures. It, therefore, ensures safe surgical intervention.

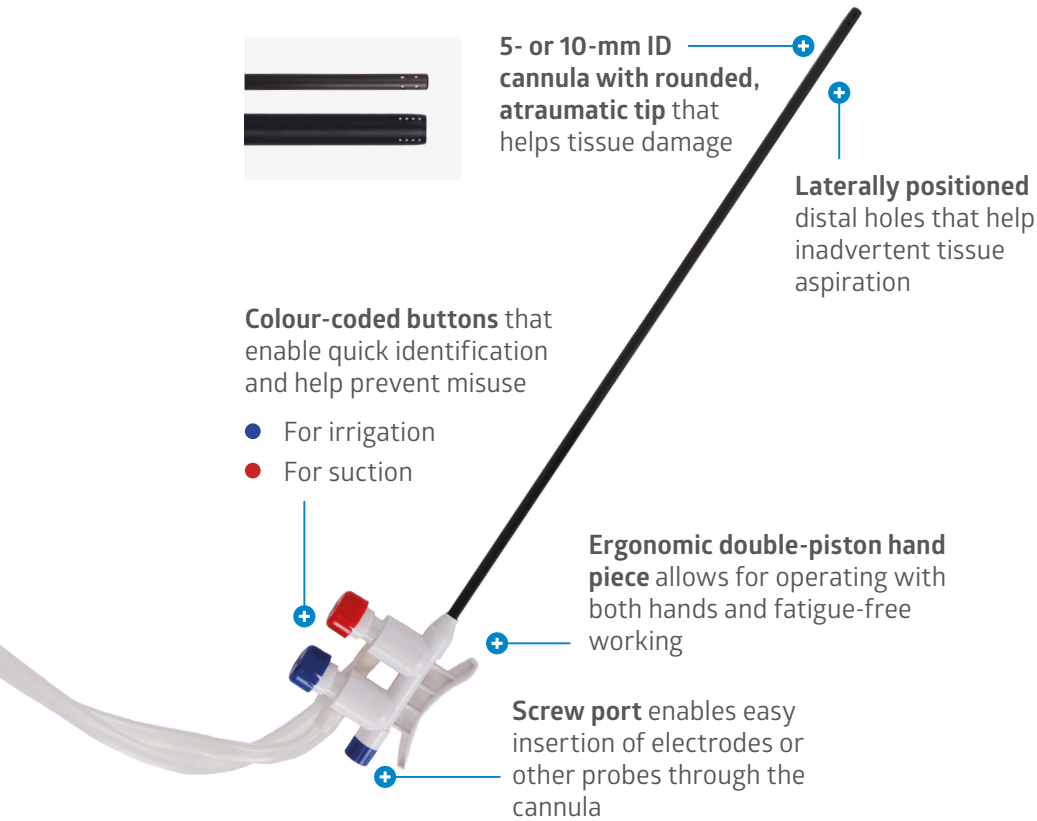

Benefits and features:

- DEHP-free
- Low torque tubing with easy-peel separation
- 33-cm cannula suction-irrigation set, or 45-cm for bariatric and robotic surgery
- 45-cm single cannula available
- Single or double spike to connect to 1 or 2 irrigation bags



Cannula made of **non-conductive polycarbonate**;^{*} compatible with electrosurgical procedures

^{*}5-mm cannula only

Item code	Description	Suction-irrigation cannula			Packaging
		Outer diameter	Length	Lateral holes	
OR1200S	Suction-irrigation set with single spike	5 mm	33 cm	8	12/cs
OR1201D	Suction-irrigation set with double spike	5 mm	33 cm	8	12/cs
OR1245S	Suction-irrigation set with single spike	5 mm	45 cm	8	12/cs
ORC45	Suction-irrigation cannula	5 mm	45 cm	8	10/cs
ONG299601CST12	Suction-irrigation set with single spike	10 mm	33 cm	8	5/cs

MED-LAP Extraction Bags

For the removal of tissue or organ specimens

Medline offers two types of extraction bags with unique features and a degree of sophistication. The Extraction Bag with Pusher is an intuitive product, which is easy to handle and widely used during laparoscopic procedures.



+ **Insertion sleeve with pusher**
To help insert the extraction bag into the abdomen



+ **Self-opening wire**
For easy insertion of the specimen

+ **Polyurethane bag**
Tear-resistant material



Benefits and features:

- Suitable for use in combination with a standard 10 or 12-mm trocar


Item code	Description	Packaging
100200	MED-LAP Extraction Bag 200 ml	5/box
100400	MED-LAP Extraction Bag 400 ml	5/box
100800	MED-LAP Extraction Bag 800 ml	5/box



MED-LAP Extraction Bag with Grip

Steerable, self-opening, detachable

The Extraction Bag with Grip is a sophisticated product that clinicians can handle easily during laparoscopic procedures.

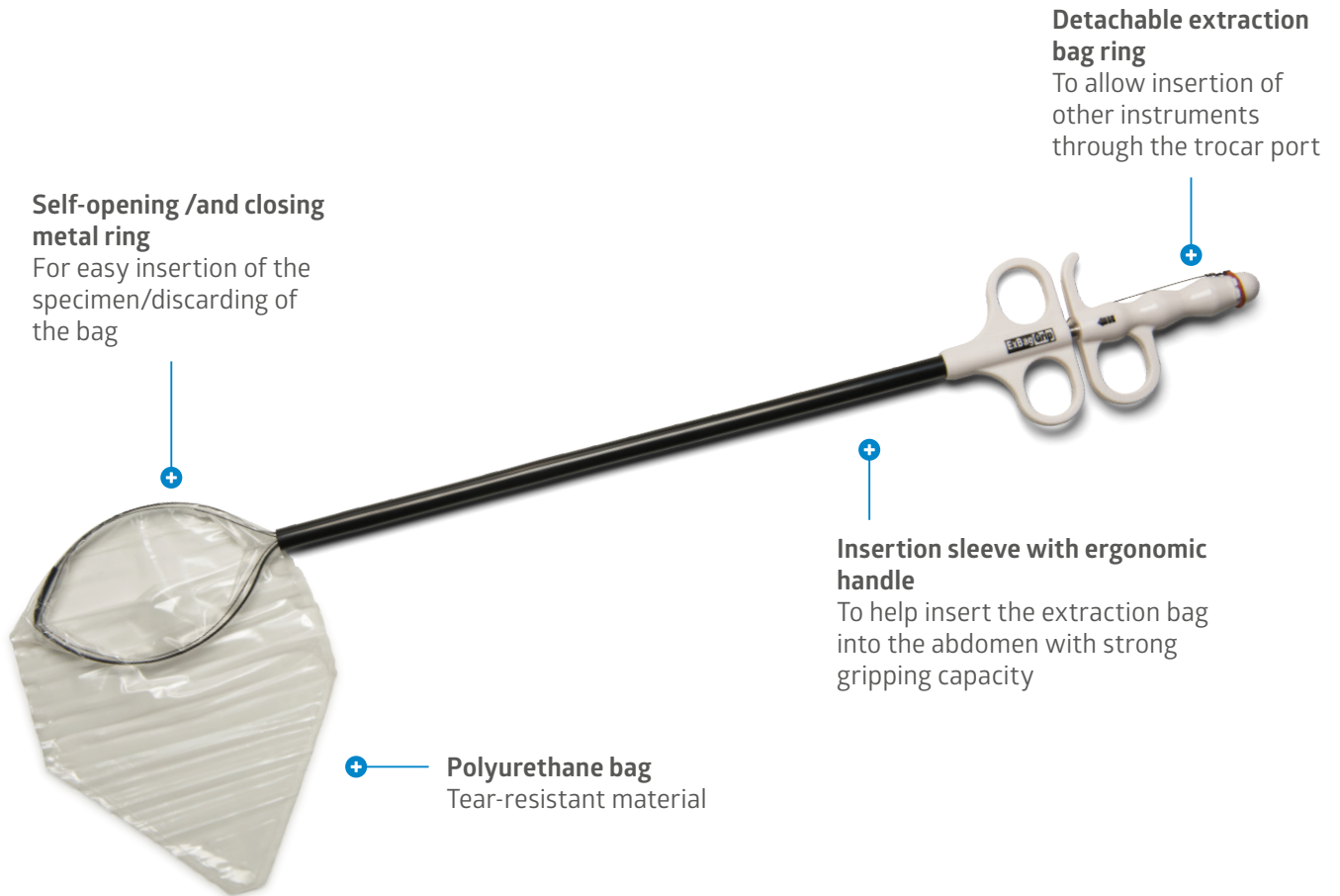

Self-opening /and closing metal ring

For easy insertion of the specimen/discarding of the bag

Detachable extraction bag ring

To allow insertion of other instruments through the trocar port

Polyurethane bag
Tear-resistant material

Insertion sleeve with ergonomic handle

To help insert the extraction bag into the abdomen with strong gripping capacity



Benefits and features:

- Suitable for use in combination with a standard 10 or 12-mm trocar

Item code	Description	Packaging
0110310	MED-LAP Extraction Bag With Grip 300 ml	5/box
0110320	MED-LAP Extraction Bag With Grip 550 ml	5/box
0110380	MED-LAP Extraction Bag With Grip 800 ml	5/box



For more information on our products, please contact your Medline account manager or visit our website: www.medline.eu/uk



Medline Industries Ltd.
3rd Floor
Quayside Wilderspool Business Park
Greenalls Avenue
Warrington WA4 6HL
United Kingdom
Tel.: +44 844 334 5237
www.medline.eu/uk
uk-customerservice@medline.com

Medline Ireland Ltd.
1 Grant's Row
Mount Street Lower
Dublin 2
Ireland
Tel.: +353 1 691 73 06

www.medline.eu/ie
ie-customerservice@medline.com



FOLLOW US  

ONGxxx set is a class IIa sterile medical device.  0459 (G-MED)  Rocamed France,
ZI de Signes, Allée de Stockholm
83870 Signes, France

ORISxxx are class IIa sterile medical devices.  2460 (DNV Product Assurance AS)  OBELIS S.A.
Bd. Général Wahis, 53
1030 Brussels, Belgium

The Veress needles are class IIa sterile medical devices.  1639 (SGS Belgium NV)  Shanghai International Holding Corp. GmbH (Europe)
Erfenstrasse 80
20537 Hamburg, Germany

The anti-fog solution is a class IIa sterile medical device.  2797 (BSI, The Netherlands)  Emergo Europe
Prinsessegracht 20
2514 AP The Hague, The Netherlands

OR12xxx and ORC45 are class IIa sterile medical devices.  0197 (TÜV Rheinland)  Luxus Lebenswelt GmbH
Kochstraße 1
47877 Willich, Germany

These products are intended to be used by healthcare professionals. Before use, consult instructions and precautions on the corresponding labelling.

We reserve the right to correct errors that may occur within this brochure.

© 2023 Medline is a registered trademark of Medline Industries, LP. Three Lakes Drive, Northfield, IL 60093, USA. ML1048_EN 08/2023