

# Sahara Absorbent Fitted Sheet

Improve efficiency. Promote safety. Reduce costs.



The new Sahara Absorbent Fitted Sheet is designed to help emergency departments cover trolleys in A&E and ambulances. The fitted design guarantees fast and easy application while the non-slip back sheet keeps the sheet in place, even during rapid patient transportation.

## Key benefits and features


- Ensures patient comfort with its quilted, smooth outer layer
- Maintains skin integrity, thanks to core polymer technology that rapidly wicks moisture away
- Impervious, non-slip back sheet
- Stays in place, thanks to the elasticated band around the edges
- Reduces the risk of cross-contamination because it is single-use
- Minimises the risk of allergic reactions because it is latex-free

### Outer layer

- Soft, non-woven; helps improve patient comfort

### Middle layer

- Super-absorbent polymer rapidly absorbs and locks in moisture
- Quilted pattern rapidly wicks moisture away from the surface

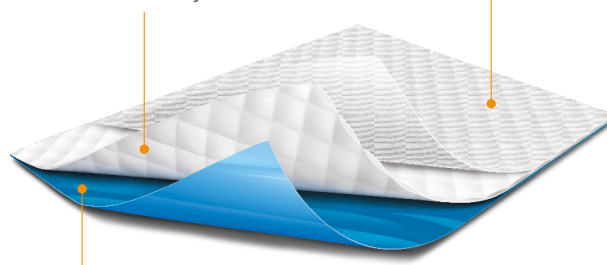

### Bottom layer

- Impervious polyethylene (PE) backing provides barrier protection; helps prevent fluids from compromising operating table pads over time

## Ordering information

Item number	Description	Quantity
DYNDE4090FTD	Sahara Fitted 228 x 101 cm	20/case

**For more information** on this product, please contact your Medline account manager or visit our website: [uk.medline.eu](http://uk.medline.eu)

Sahara Absorbent Fitted Sheet is a class I NS medical device intended to be used by healthcare professionals. Before use, consult instructions and precautions on the corresponding labelling. We reserve the right to correct errors that may occur within this brochure.

© 2020 Medline Industries, Inc. Medline is a registered trademark of Medline Industries, Inc. ML1070\_EN 11/2020.

Medline International France SAS  
5, rue Charles Lindbergh  
44110 Châteaubriant, France

