

Signature™ Latex Essential

Sterile latex powder-free surgical gloves



- Made from natural rubber latex
 - » **Caution:** Safe use of these gloves by or on latex-sensitised individuals has not been established
 - » **Caution:** This product contains natural rubber latex, which may cause allergic reactions, including anaphylactic response
- Powder-free to help reduce powder-related skin irritation and post-operative complications associated with granulomas
- Exclusive Medline SH8 glove mould with an independent thumb design that allows for an anatomical fit and more natural movement in the fingers, thumb and palm
- Thinner than Signature™ Latex Surgical Gloves for enhanced tactile sensitivity
- Short-dip technology does not apply to E-Z Glide polymer coating near the cuff, which helps reduce cuff roll down

Specifications

Thickness:



Finish:



Texture:

Smooth finish for enhanced tactile sensitivity

Colour:



Benefits and features



Ideal for general surgery



Effortless donning with damp or dry hands is possible with included E-Z Glide polymer coating



Cuff roll-down (common issue with many gloves) is reduced, thanks to anti-roll-down cuff (ARC) technology



Recyclable and waste-reducing packaging



Enhanced comfort through improved design



This product is made with natural rubber latex



Excellent choice as a top glove when double-gloving



Ordering information

Item code	Size
MSG5955	5 ½
MSG5960	6
MSG5965	6 ½
MSG5970	7
MSG5975	7 ½
MSG5980	8
MSG5985	8 ½
MSG5990	9


50 pairs/box, 200 pairs/case

Length and width


Item code	Length	Width
MSG5955	276 mm	69 mm
MSG5960	280 mm	78 mm
MSG5965	285 mm	79 mm
MSG5970	285 mm	89 mm
MSG5975	285 mm	93 mm
MSG5980	286 mm	101 mm
MSG5985	288 mm	107 mm
MSG5990	287 mm	108 mm

E-Z Glide

E-Z Glide polymer coating uses triple-layered technology, including a proprietary surfactant coating, to ensure easy donning with damp or dry hands.

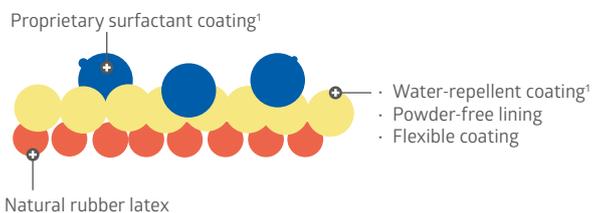

Regulations and quality standards

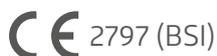

Medline International France SAS
5 Rue Charles Lindbergh
44110 Châteaubriant, France

- Manufacturer EN ISO 13485-compliant
- Class IIa medical device; compliant with Medical Device Directive 93/42/EEC
- Category III personal protective equipment (PPE); compliant with PPE Regulation (EU) 2016/425

Learn more. Contact your Medline account manager or visit: www.medline.eu/uk

This glove is a class IIa sterile medical device intended to be used by healthcare professionals. Before use, consult instructions and precautions on the corresponding labelling.

All rights reserved. We reserve the right to correct errors that may have occurred.

© 2023 Medline Industries, LP. Medline is a registered trademark of Medline Industries, LP. Three Lakes Drive, Northfield, IL 60093, USA. ML572_EN 07/2023