




CENTURION
MEDICAL PRODUCTS

SorbaView® 2000


Window Dressing

**Designed to go
on easier and
last longer**

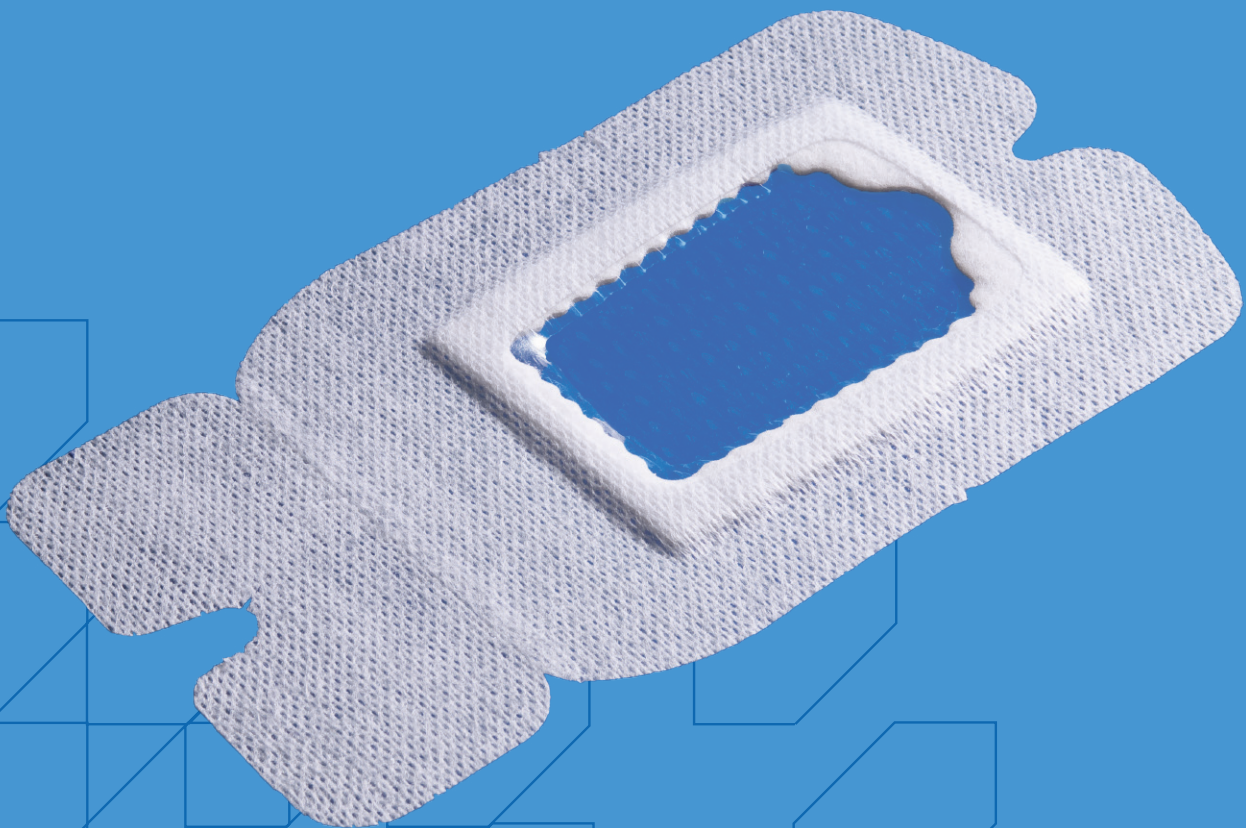

Benefits for caregivers, hospitals and patients


- Integrated pad wicks away moisture, keeping site dry
- Stays securely in place for long wear times that meet hospital protocols (up to 7 days for central lines, and 96 hours or until clinically indicated for peripheral catheters)
- Conforms and responds to body contours
- Easy to apply and remove
- Exceptionally high moisture vapour transmission rate (MVTR)
- Non-irritating, non-sensitising and non-cytotoxic adhesive
- Continuous site visibility thanks to its transparent window
- Multilayer construction that is protective, breathable and comfortable to wear
- Creates an occlusive barrier against bacteria, viruses and other contaminants, but not to moisture vapour or other normal gases
- Does not contain natural rubber latex

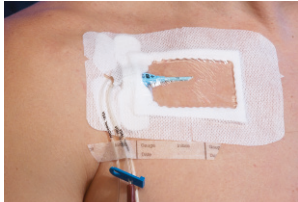

Subclavian (SV40XT)



Peripheral (SV30XT)

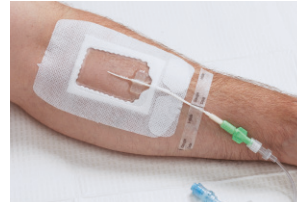


PICC (SV48VXT)



Implanted port (SV55XT)

SorbaView 2000 sizes available

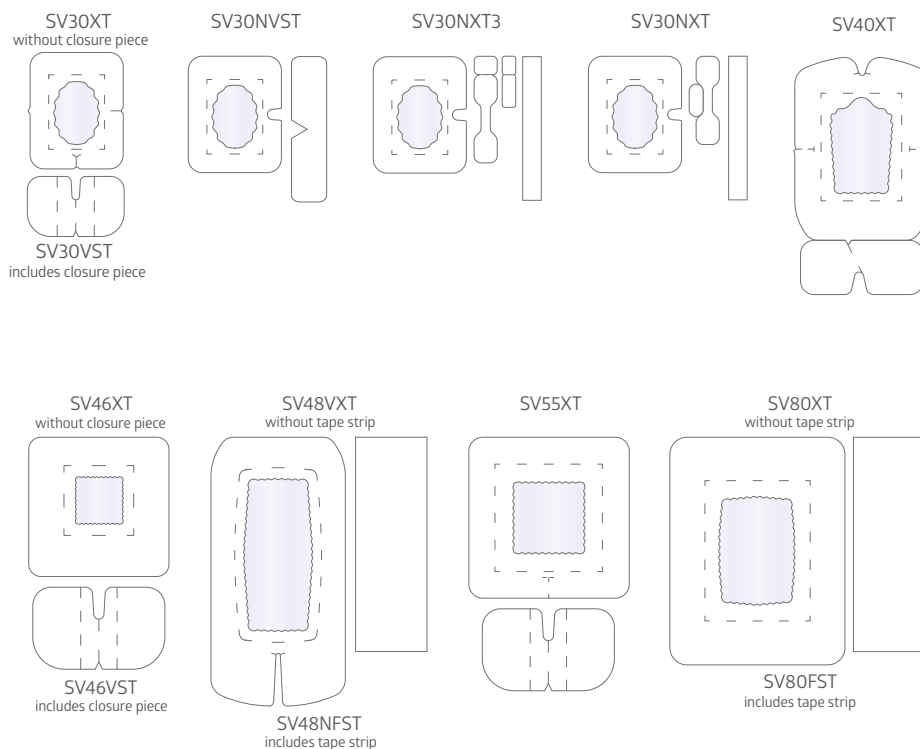

SorbaView 2000 Window Dressing sets the standard of excellence in catheter site management and securement.

SorbaView is easy to apply, and remains intact for up to seven days on central lines, reducing dressing changes, and thus lowering the risk of infection.^{1,2}

SorbaView allows skin to breathe while providing a proven barrier to external liquids and bacteria. Its signature design allows vapours to escape through the film window while the absorbent pad wicks away excess moisture. This prevents build-up of fluid at the site and keeps moisture from breaking down the adhesive under the fabric border. The adhesive stays dry and dressings stay on!



Medline Industries Ltd
3rd Floor
Quayside Wilderspool Business Park
Greenalls Avenue
Warrington WA4 6HL
United Kingdom
Tel.: +44 844 334 5237
www.medline.eu/uk
uk-customerservice@medline.com

Medline Ireland Ltd
1 Grant's Row
Mount Street Lower
Dublin 2
Ireland
Tel.: +353 1 691 73 06
www.medline.eu/ie
ie-customerservice@medline.com

FOLLOW US  


A comprehensive line of Centurion catheter securement products is available to further help with stabilising and protecting wound sites.

¹Lee, M. Seven Day Dwell Times for Vascular Access Device Exit Site Dressings, Lee Medical Inc. 2005. ²Garcia, R et al. Three Years Experience in Implementing HICPAC Recommendations for the Reduction of Central Venous Catheter-Related Bloodstream Infections. Brookdale University Medical Center (BUMC), 2003.

SorbaView 2000 dressings are class I sterile medical devices intended to be used by healthcare professionals. Before use, consult instructions and precautions on the corresponding labelling.



Medline International France SAS
5 Rue Charles Lindbergh
44110 Châteaubriant, France

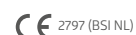

All rights reserved. We reserve the right to correct errors that may have occurred.

© 2022 Medline Industries, LP. Medline is a registered trademark of Medline Industries, LP. ML943_EN 06/2022

© 2022 Centurion Medical Products Corporation LIT063V3E