




Secure protection everywhere you need it

Suresite® ensures a fixed, comfortable fit at every application site.

- Short peripheral and midline venous catheters
- Central venous catheters (CVCs)
- Peripherally inserted central catheters (PICCs)
- Tunnelled CVCs
- Implanted infusion devices
- Arterial catheters
- Dialysis catheters

Ordering information

Item number	Description	Packaging
VAD1680-INT	4.4 x 3.8 cm	400/cs; 100/bx
VAD1682-INT	5.7 x 5.1 cm	400/cs; 100/bx
VAD1882-INT	7 x 8.9 cm	200/cs; 50/bx
VAD1683-INT	7 x 7.3 cm	400/cs; 100/bx
VAD1685-INT	8.9 x 11.4 cm	200/cs; 50/bx
VAD1688-INT	10.2 x 11.4 cm	200/cs; 50/bx
VAD1689-INT	10.2 x 15.5 cm	100/cs; 25/bx

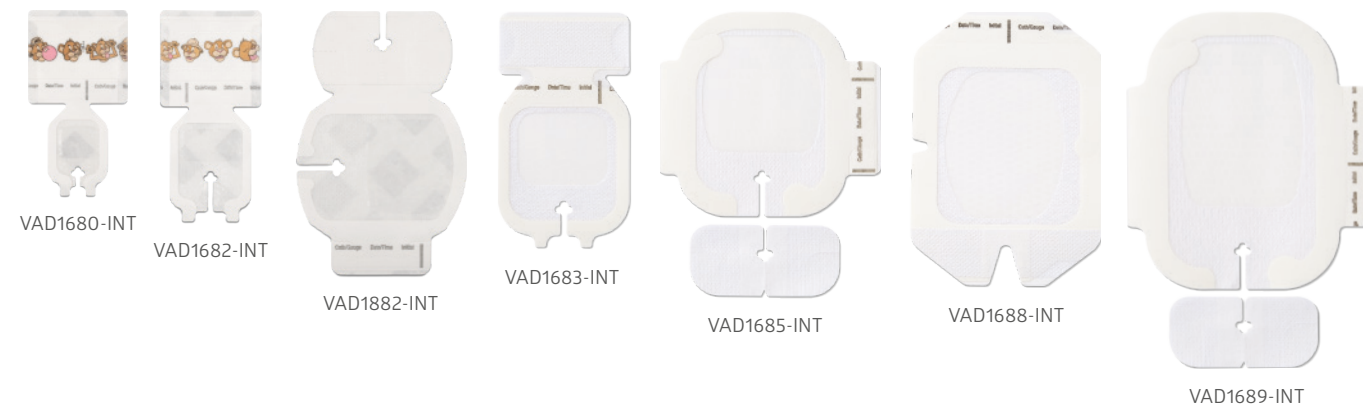

Hold it steady. Keep it clean.

Keep catheters sterile and in place.

Suresite® Select Sterile Dressings improve the patient's comfort and simplify the caregiver's work. The dressings form a complete seal, hold the catheter firmly in place and help prevent harmful microorganisms from entering the site. Suresite's exceptionally wide viewing window also makes it easy to monitor for signs of infection.

Ensure a sure fit. Contact your Medline account manager for information on adding Suresite to a custom kit.



Medline Industries Ltd
3rd Floor
Quayside Wilderspool Business Park
Greenalls Avenue
Warrington WA4 6HL
United Kingdom
Tel.: +44 844 334 5237
Fax: +44 844 334 5238
www.medline.eu/uk
uk-customerservice@medline.com

FOLLOW US



Protection you can count on

Suresite® Select Dressings are class I sterile medical devices intended to be used by healthcare professionals. Before use, consult instructions and precautions on the corresponding labelling.



We reserve the right to correct errors that may occur in this brochure.
© 2021 Medline is a registered trademark of Medline Industries, Inc. Three Lakes Drive, Northfield, IL 60093, USA. ML1045_EN 06/2021.

Transparent film
provides a breathable,
sterile barrier to
external contaminants

Innovative adhesives
ensure a secure hold
and manage moisture

Frame delivery
ensures accurate,
easy placement
and gentle
removal

Wide viewing window
makes it easy to monitor
for signs of infection

Clover notch fits easily
around and conforms to
hubs and lumens

How to apply

1. Remove liner from dressing to expose adhesive
2. Position notch edge of dressing over catheter hub
3. Peel frame away from border edge while applying pressure to border edge
4. Apply securement tape strip's pre-cut notch underneath catheter's moulded junction



VAD1683-INT 7 x 7.3 cm



VAD1685-INT 8.9 x 11.4 cm



VAD1688-INT 10.2 x 11.4 cm



VAD1689-INT 10.2 x 15.5 cm




How to remove

1. Remove all securement tape strips from dressing
2. Gently grasp soft cloth edges below keyhole notch
3. Lift dressing on each side from the skin using low-and-slow technique
4. Peel dressing towards catheter insertion site until removed



VAD1683-INT 7 x 7.3 cm

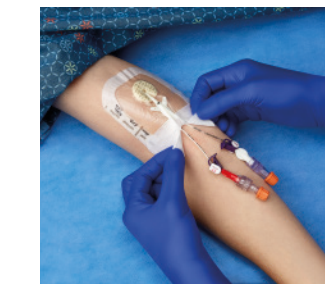

VAD1685-INT 8.9 x 11.4 cm



VAD1688-INT 10.2 x 11.4 cm



VAD1689-INT 10.2 x 15.5 cm

