

Suction Cannulas and Yankauers

A complete selection of suction system accessories



**Trusted fluid
management**

Tools designed for every need

Medline's suction instruments come in a variety of styles to suit your preferences and clinical needs for the removal of fluids during surgical procedures. They deliver performance and manoeuvrability to provide optimal control and visibility while helping to minimise the risk of tissue trauma. All suction instruments are latex-free and have been double-wrapped individually to ensure sterility. Depending on the need, Medline offers flexible and rigid suction instruments.

Flexible yankauer with female connector

These yankauers with a female connector offer an efficient suctioning solution. Their flexibility facilitates suctioning in deep and sensitive areas without damaging tissues.

Each yankauer is made of PVC and is suitable for use during oral aspiration, suction-assisted laryngoscopy and airway decontamination.

Code	Description	Length	Diameter	Pack factor
Y1800F	Flexible yankauer with female connector	250 mm	18 CH	50
Y2200F	Flexible yankauer with female connector	250 mm	22 CH	50

Rigid classic yankauer

Medline offers a full range of common, rigid and single-use suction cannulas for suction of fluids, with a wide range of open tips (fine, regular, large) to meet your requirements for all applications. The range includes lateral eyes for controlled atraumatic suction, with a high-quality finish and heat-moulded rounded distal tips and bevelled edges to minimise tissue trauma, even on the most delicate of tissues.

Code	Description	Length	Diameter	Holes	Vacuum control	Pack factor
Y1200	Fine tip; for paediatric, dental and ENT use	270 mm	12 CH	2	No	50
Y1201		270 mm	12 CH	2	Yes	50
Y1500		200 mm	15 CH	0	No	100
Y1501	Regular tip; for pharyngeal use	200 mm	15 CH	0	Yes	100
Y1800		250 mm	18 CH	6	No	50
Y1801		250 mm	18 CH	6	Yes	50
Y2200	Large tip; for orthopaedic use, large particles and significant volumes of fluid	270 mm	22 CH	6	No	50
Y2201		270 mm	22 CH	6	Yes	50
Y2800		270 mm	28 CH	4	No	50
Y2801		270 mm	28 CH	4	Yes	50

Rigid yankauer with tapered bulb tip and ergonomic handle

Code	Description	Length	Diameter	Holes	Vacuum control	Pack factor
Y80	Tapered bulb tip	262 mm	12 CH	0	No	50
Y82	Tapered bulb tip	262 mm	12 CH	0	Yes	50

Rigid yankauer with flexible connector

This yankauer with a flexible connector is manufactured from transparent PVC and can be used for general suction purposes, dental suction and orthopaedic procedures. The flexible graduated proximal end ensures a universal tubing connection.

Code	Description	Length	Diameter	Holes	Vacuum control	Pack factor
P2000	Straight tip	280 mm	20 CH	4	No	50
P20FT	Flat tip/duckbill	280 mm	20 CH	4	No	50

Rigid yankauer with ergonomic handle

Shatter-resistant PVC can be used for moderate organ retraction or tissue dissection. The ergonomic, non-slip grip handle ensures comfortable handling during long procedures.

Code	Description	Length	Diameter	Holes	Vacuum control	Pack factor
Y2400	Straight tip	253 mm	24 CH	4	No	50
Y240B	Bulb tip	260 mm	24 CH	4	No	50

Rigid yankauer with ON/OFF switch handle

The switch system allows the operator to manually open and close the vacuum to control the suction level and reduce noise. The anatomically designed handle makes usage more comfortable. This yankauer also comes with a graduated proximal end for a secure universal tubing connection.

Code	Description	Length	Diameter	Holes	ON/OFF vacuum control	Pack factor
Y240FS	Straight tip	297 mm	24 CH	4	Yes	50
Y240FB	Bulb tip	297 mm	24 CH	4	Yes	50

Poole suction cannula

This cannula is made of a two-piece fixed construction:

- The inner cannula has 4 side holes and an open distal tip to promote smooth flow.
- The outer flexible sheath has multiple side holes and an open distal end for dispersed suction, without the risk of tissue invagination.

The Poole cannula is a specialty suction instrument, designed to evacuate pooled blood, fluid and debris during deep abdominal surgeries, including caesarean sections.




Code	Description	Length	Diameter	Holes	Vacuum control	Pack factor
Y20PS	Flexible straight tip	297 mm	20 CH	20 + 4	No	50



Yankauer sets


These sterile yankauer suction handles are pre-attached to 1.8- or 3-metre suction tubing with ribbed female connectors.

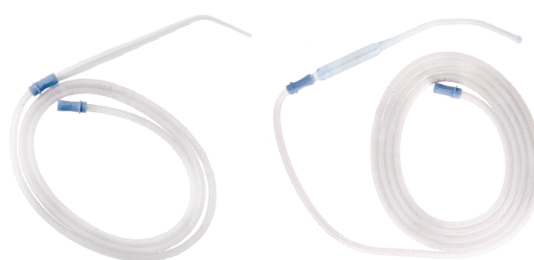

Code	Description	Yankauer				Suction tubing		Pack factor
		Diameter	Length	Tip holes	Vacuum control	Internal diameter	Length	
ORNEX66AY1200	Fine tip yankauer set	12 CH	270 mm	2	No	6 mm	1.8 m	30
ORNEX610Y1800	Regular tip yankauer set	18 CH	250 mm	6	No	6 mm	3 m	30
ORNEX610Y2200	Large tip yankauer set	22 CH	270 mm	6	No	6 mm	3 m	30
ORNEX610Y80	Fine and tapered bulb tip yankauer set; with ergonomic handle	12 CH	280 mm	0	No	6 mm	3 m	30

Orthopaedic suction set

The orthopaedic suction set with a filter helps maintain a clear operative field. The special filter built into the handle automatically traps solids and can be exchanged in a matter of seconds when filled. Increased filter capacity makes the suction set ideal for all procedures that produce a large amount of debris.

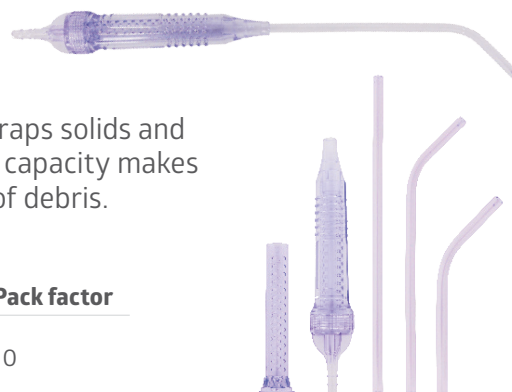

Code	Description	Pack factor
YORTH0	Orthopaedic suction set; with 1 ribbed handle, 1 curved tip (7 x 230 mm), handle with 2 filters and 3 possible tips	10





Medline Industries Ltd
3rd Floor
Quayside Wilderspool Business Park
Greenalls Avenue
Warrington WA4 6HL
United Kingdom
Tel.: +44 844 334 5237
www.medline.eu/uk
uk-customerservice@medline.com

Medline Ireland Ltd
1 Grant's Row
Mount Street Lower
Dublin 2
Ireland
Tel.: +353 1 691 73 06

www.medline.eu/ie
ie-customerservice@medline.com

FOLLOW US  

All yankauers and yankauer sets are class IIa sterile medical devices intended to be used by healthcare professionals. Before use, consult instructions and precautions on the corresponding labelling.

References YxxxxF:  Shanghai International Holding Corp. GmbH (Europe)  0197 (TÜV Rheinland LGA Products GmbH)

We reserve the right to correct errors that may occur in this brochure.

© 2023 Medline is a registered trademark of Medline Industries, LP, Three Lakes Drive, Northfield, IL 60093, USA. ML468_EN 01/2023.



Luxus Lebenswelt GmbH
Kochstr. 1
47877 Willich, Germany



0197 (TÜV Rheinland LGA Products GmbH)