



FLUID MANAGEMENT



MED-RIGID

Disposable Rigid Suction Canisters

MED-RIGID Suction Canisters

Discover strong, durable construction and ease of use.

Disposable suction system based on a canister/holder combination for the collection of bodily fluids, liquids and aspiration waste. The canister holder is reusable and the rigid canister is single patient use. The MED-RIGID system features single-use disposable hard locking lid attachments and a clear transparent canister with ability to stand alone. This device encourages the proper disposal of infectious liquid medical waste and enhances worker safety.

Special brackets and stands are available for convenient mounting configurations and mobility.

The patient port is equipped with a 90° elbow connector. It can be removed or used to facilitate the connection and avoid kinking of the patient tubing.

Durable locking dome lid prevents leaking and spilling

>99.99% rated Bacterial Filtration Efficient (BFE) filter to avoid vacuum source contamination

Hard and flat stand alone bottom




Transparent, shatter-resistant and fully disposable

Accurate graduation markings

MED-RIGID Canister Holders

Reusable Canister Holders are used to support the MED-RIGID canister in any slide mount such as the wall support or mobile trolley. Built with preattached red vacuum tube for connection to MED-RIGID canister and ON/OFF switch for vacuum control.

Item No.	Description	Pkg.
OR159	MED-RIGID Canister Holder, for 1.2L canister	20/cs
OR139	MED-RIGID Canister Holder, for 2 or 3L canister	20/cs

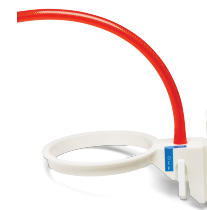

MED-RIGID Canisters

MED-RIGID Canisters are single patient use. Available in 1.2, 2L, 3L and even 12L each featuring a hard, disposable body and lid.

Item No.	Description	Pkg.
OR212	MED-RIGID Canister, 1.2L	40/cs
OR220	MED-RIGID Canister, 2L	40/cs
OR230	MED-RIGID Canister, 3L	40/cs
OR12L	MED-RIGID Canister, 12L For this specific canister, please use its dedicated metallic trolley (Item No. OR12LTR)	4/cs



MED-RIGID Canisters with Suction Tubing

MED-RIGID Canisters are available in kits, including a non-sterile ORNEX tubing. ORNEX is a low memory flexible tubing with ribbed external surface to increase its resistance.

Item No.	Description	Pkg.
OR312	MED-RIGID Canister 1.2L + 5mm x 1.8m Tubing	30/cs
OR320	MED-RIGID Canister 2L + 5mm x 1.8m Tubing	30/cs
OR319	MED-RIGID Canister 3L + 6mm x 3.0m Tubing	30/cs



MED-RIGID Tandem Tubing

Used for the connection of several canisters and automatic transfer of fluids from one canister to the next one.

Item No.	Description	Pkg.
OR542	Tandem Tubing for MED-RIGID Canister	50/cs



HARDWARE and Dedicated Equipment

Disposable Specimen Traps


The disposable specimen traps are used to collect and retain samples for analysis. They can be fixed on the accessory port of the MEDLINE canisters and connected to a 7mm or 10mm suction tubing.

Item No.	Description	Pkg.
OR123	Short specimen sock (9cm), with male connector diameter 7mm	100/cs
OR125	Standard specimen sock (23cm), with male connector diameter 7mm	100/cs
OR124	Large specimen sock (23cm), with male connector diameter 10mm	100/cs
CA6000	Specimen cup, diameter 10cm	10/cs



Mounting Brackets

Item No.	Description	Pkg.
DYNDSCW	Plastic Wall Mount	20/cs
LCLC1426CE	Metal Rail Clamp	1ea/cs

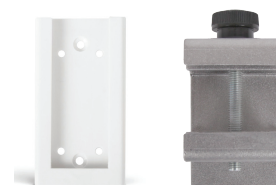

Trolleys

Enable advanced and customised configurations.

Item No.	Description	Pkg.
DYNDSCC	Metallic Canister Trolley Adjustable height between 61 and 106cm Can hold 4 canisters	1/cs
ORBASE	Canister Trolley Fixed height 50cm. Can hold 2 canisters	1/cs



Manifold Tubing

Our reusable manifold tubings enable customised set-up by connecting a multiple canister network to the vacuum source.

Item No.	Description	Pkg.
DYND2MNF	Manifold Tubing, for 2 canisters	10/cs
DYND3MNF	Manifold Tubing, for 3 canisters	10/cs
DYND4MNF	Manifold Tubing, for 4 canisters	10/cs



Medline Industries Ltd
No 5 Booths Park
Chelford Road
Knutsford, Cheshire WA16 8GS
England
Tel: +44 844 334 5237
Fax: +44 844 334 5238
www.medline.eu/uk
uk-customerservice@medline.com